

DIRECTA
FenderWedge®

DIRECTA
FenderMate®

Protect.
Prepare.
Restore.

EN

DIRECTA FenderWedge®

A wedge and steel shield in one piece, for use between adjacent teeth during preparation

Accidental bur contact with a neighbouring tooth is a familiar risk in everyday practice. FenderWedge pairs a steel plate with a plastic wedge in a single piece, designed to help shield the adjacent tooth and gingiva during preparation. The same combination is intended to ease matrix placement afterwards, so the workflow can continue smoothly from preparation to restoration.



Shield for the Adjacent Tooth

The steel plate acts as a barrier between bur and neighbouring tooth, and is designed to help reduce the risk of accidental damage.

Gentle on Gingiva


Shield and wedge work together and are intended to help reduce trauma to the gingival tissue while you prepare.

One Placement, Several Roles

Placed once, the same unit supports your preparation and remains in position for the next step.

Ready for the Matrix

The wedge is designed to provide pre-separation, helping to facilitate matrix placement and proximal contact formation.

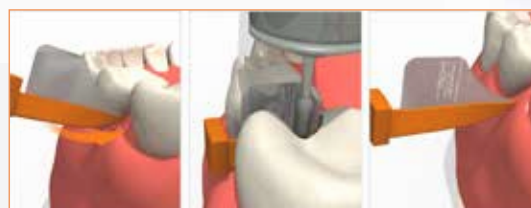

Prod No Product Description

602750	FenderWedge Assortment, 4x36 pcs
602751	FenderWedge Small, orange, 36 pcs
602752	FenderWedge Medium, green, 36 pcs
602753	FenderWedge Large, yellow, 36 pcs
602754	FenderWedge X-small, lilac, 36 pcs
602800	FenderWedge Value Pack, orange, 100 pcs
602801	FenderWedge Value Pack, green, 100 pcs
602802	FenderWedge Value Pack, yellow, 100 pcs
602803	FenderWedge Value Pack, lilac, 100 pcs
602679	FenderWedge & FenderMate Intro-Kit*, 8 & 20 pcs

*The Intro-Kit contains 3 FenderWedge X-small, 3 Small, 1 Medium and 1 Large, together with 8 FenderMate Left Narrow, 2 Left Regular, 8 Right Narrow and 2 Right Regular.



DIRECTA

FenderWedge® Prep

Designed to help protect the adjacent tooth and papilla in crown, veneer and inlay preparations

FenderWedge Prep is a plastic wedge with an integrated, perforated steel plate, made for larger preparations such as crowns, veneers and inlays. Its thin, boat-shaped tip is intended to enter tight proximal spaces and gently compress the gingiva, which may improve access while helping to limit tissue trauma. The plate can be bent before insertion to rest against the adjacent tooth, supporting access and visibility.



Adjacent Tooth Shield

The steel plate is designed to help reduce the risk of accidental bur damage to the neighbouring tooth.



Boat-shaped Tip

The slim tip is designed to slip into tight interproximal spaces and compress the gingiva, supporting access during subgingival preparation.

Bend to Fit

Shape the plate before insertion so it rests against the adjacent tooth, helping to support access and visibility.

Papilla Care

Designed to help protect the papilla during crown and veneer work.

Prod No	Product Description
602749	FenderWedge Prep, grey, 36 pcs
602811	FenderWedge Prep Value Pack, grey, 100 pcs



DIRECTA FenderMate®

Sectional matrix for Class II composite restorations

FenderMate is a pre-curved sectional matrix and wedge in one piece, designed to support Class II composite restorations, proximal contact formation and cervical adaptation. It can be placed from either the buccal or lingual side. The matrix runs from the base of the wedge to just above the occlusal surface, and an angled wing on the wedge faces the adjacent tooth.



Contact and Contour

The pre-shaped contact is designed to adapt to the proximal surface, supporting contact formation and anatomical contouring.

One-piece design

One item to place instead of two, which is intended to simplify setup.

Wing for Separation

The wing applies pressure that is intended to help maintain separation and support adaptation at the cervical margin, which may help reduce the risk of overhangs.

Interproximal Placement

Designed for insertion into the interproximal space with minimal preparation of the site.

Prod No Product Description

602770	FenderMate Assortment, 4x18 pcs
602771	FenderMate left regular, green, 18 pcs
602772	FenderMate left narrow, light green, 18 pcs
602773	FenderMate right regular, blue, 18 pcs
602774	FenderMate right narrow, light blue, 18 pcs
602804	FenderMate Value Pack, left regular, green, 100 pcs
602805	FenderMate Value Pack, left narrow, light green, 100 pcs
602806	FenderMate Value Pack, right regular, blue, 100 pcs
602807	FenderMate Value Pack, right narrow, light blue, 100 pcs
602679	FenderWedge & FenderMate Intro-Kit*, 8 & 20 pcs


*The Intro-Kit contains 3 FenderWedge X-small, 3 Small, 1 Medium and 1 Large, together with 8 FenderMate Left Narrow, 2 Left Regular, 8 Right Narrow and 2 Right Regular.



DIRECTA **FenderMate® Temp**

Matrix and protection for temporary restorations

FenderMate Temp combines a stainless-steel matrix and a plastic wedge for preparing and filling temporary proximal restorations. It can be placed in the proximal area before preparation and is designed to help reduce the risk of accidental damage to the adjacent tooth. It is also intended to help limit blood from the gingival pocket in the treatment area. After preparation, the matrix can be curved to follow the tooth anatomy, supporting precise placement of the temporary filling and helping to reduce the risk of overhangs.

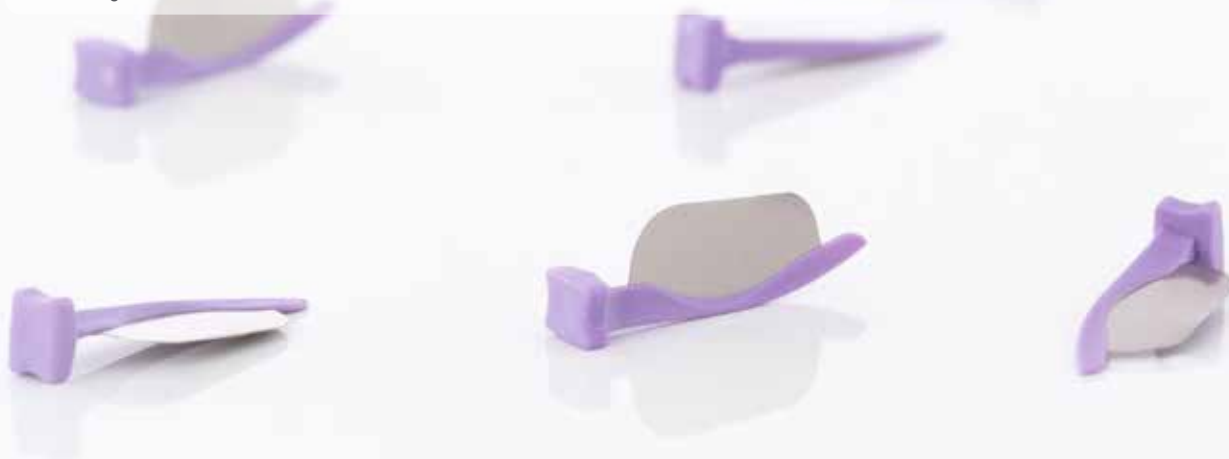

Made for Temporary Fillings

Intended for temporary proximal restorations, from preparation to filling.

Placed Before Preparation


Inserted in advance, it is designed to help shield the adjacent tooth and limit blood from the gingival pocket in the working area.

Curve to the Tooth

The matrix can be shaped around the tooth to follow its anatomy, which may help minimise the need for finishing.

One-piece Design

Two parts in one: one item serves as both protection during preparation and matrix during filling.

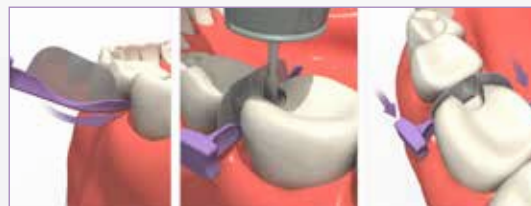


Prod No	Product Description
602786	FenderMate Temp, lilac, 18 pcs
602813	FenderMate Temp Value Pack, lilac, 100 pcs



DIRECTA **FenderMate® Prime**

Protective matrix for primary teeth

FenderMate Prime is a preparation shield and matrix in one, made for Class II restorations in primary teeth. A plastic interproximal wedge carries an integrated steel plate, and its boat-shaped tip is designed to compress the gingiva while helping to limit tissue disruption during placement. Two lengths are available. The long version (neon green) is bendable and can be shaped to the contour of the tooth. The short version (neon yellow) has a more rigid plate for very tight interproximal spaces.

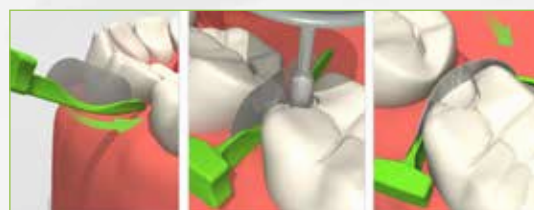

Proportioned for Primary Teeth

Sized to suit the smaller anatomy and narrower interproximal spaces of primary teeth.

Choose Your Length

Long (neon green): bendable plate for general use. Short (neon yellow): rigid plate that should not be bent, for very tight interproximal spaces.

Shape to the Tooth

The long matrix can be adapted to the curvature of the tooth, supporting an anatomical shape and helping to minimise the need for finishing.

Protect, then Restore

One piece designed to help protect the adjacent tooth during preparation and then serve as the matrix for the Class II restoration.

Prod No Product Description

602815	FenderMate Prime Assortment, 2x18 pcs
602816	FenderMate Prime Short, neon yellow, 18 pcs
602817	FenderMate Prime Long, neon green, 18 pcs
602808	FenderMate Prime Value Pack, Short, neon yellow, 100 pcs
602809	FenderMate Prime Value Pack, Long, neon green, 100 pcs



DIRECTA

FenderMate® D-Coat

Dark, matt, non-stick diamond coating

FenderMate D-Coat is a no-ring matrix system with a redesigned multilayer diamond coating, designed to reduce composite adhesion and support visual contrast with a dark, matt surface. The pre-curved matrix has a contact point and a specially designed wedge, which supports cervical adaptation and helps reduce the risk of overhangs.

Colour coding: green and blue show a left or right curve, light shades a thin wedge and dark shades a thick wedge.



Non-stick Coating

The amorphous diamond coating is designed to reduce adhesion of composite to the surface, which is intended to ease matrix removal after curing.

Dark, Matt Finish

The dark surface contrasts with composite resin, and the matt finish is designed to reduce light reflection, supporting visual contrast.

Wedge and Matrix Together

The matrix is attached to the wedge, giving a single unit to place.

Cervical Adaptation

The wing is designed to help maintain proximal separation and support adaptation at the cervical margin, which may help reduce the risk of overhangs.

Prod No	Product Description
602994	FenderMate D-Coat Assortment, 4x18 pcs
602990	FenderMate D-Coat, left regular, green, 18 pcs
602991	FenderMate D-Coat, left narrow, light green, 18 pcs
602992	FenderMate D-Coat, right regular, blue, 18 pcs
602993	FenderMate D-Coat, right narrow, light blue, 18 pcs



Design By Dentists

Directa's products are developed, tested and evaluated by our selected team of highly-qualified and renowned Swedish and international dental professionals. The aim is to make life easier for all dental professionals in their daily clinical work.

Protect. Prepare. Restore.

DIRECTA



DIRECTA AB

Finvids väg 8
194 47 Upplands Väsby, SWEDEN
+46 8 506 505 75
info@directadental.com
www.directadental.com

DIRECTA INC

64 Barnabas Road, Unit 3
Newtown, CT. 06470, USA
+1 203.491.2273
infousa@directadental.com
www.directausa.com

DIRECTA APAC

888, Moo7, Nongsarai, Pakchong,
Nakhonratchasima THAILAND 30130
infoapac@directadental.com
www.directapac.com



[directadental](https://www.facebook.com/directadental)



[directadental](https://www.instagram.com/directadental)



[DirectaAB](https://www.youtube.com/DirectaAB)

The content of this brochure, including layout and images, is protected by copyright. ©2026 Directa AB. All rights reserved.