

**FenderWedge®, FenderMate® Prime and FenderMate® D-Coat**  
**Solutions for simple, safe and fast preparation and Class II restorations.**

**FenderMate® D-Coat**

**FenderMate® D-Coat is a unique no-ring matrix system**

The pre-curved matrix includes a contact point and comes attached to a specially designed wedge which allows tight cervical sealings with no overhang. It is designed with a coating combining a better visual contrast with anti-stick properties. For small and medium restorations. FenderMate D-Coat is available in refill and assortments. The four colours allow for an easy distinction between left or right curve (green/blue), and wedge size (light = thin, dark = thick).

**Non-stick coating**

The amorphous diamond coating prevents adhesives from sticking to the surface. Pulling the matrix out after light curing is an effortless task.

**Dark & matt surface**

The dark and matt coating contrasts with the color of the composite resin and reduces light reflection. Clinicians can therefore perform restorations under excellent visual conditions.

**One piece design**

Matrix and wedge combined for easier and safer placement.

**No cervical overhang**

The wing design helps maintain proximal separation and ensures superior cervical adaptations.



Prod No	Product Description
602994	FenderMate D-Coat Assortment, 4x18 pcs
602990	FenderMate D-Coat, left regular, green, 18 pcs
602991	FenderMate D-Coat, left narrow, light green, 18 pcs
602992	FenderMate D-Coat, right regular, blue, 18 pcs
602993	FenderMate D-Coat, right narrow, light blue, 18 pcs




**Video: FenderMate D-Coat**  
Say goodbye to matrix rings!



# DIRECTA FenderWedge®

## Protects tooth and gingiva during preparations

For Class II preparations. FenderWedge is a combination of steel plate and a plastic wedge, which facilitates easy application of a matrix.

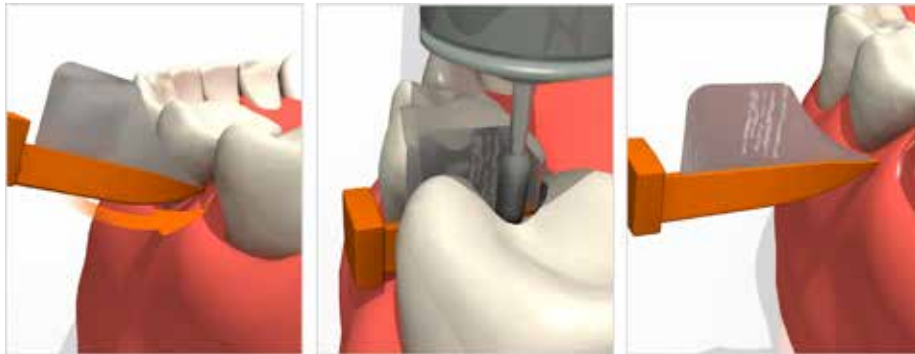

Prod No	Product Description
602750	FenderWedge Assortment, 4x36 pcs
602751	FenderWedge Small, orange, 36 pcs
602752	FenderWedge Medium, green, 36 pcs

Prod No	Product Description
602753	FenderWedge Large, yellow, 36 pcs
602754	FenderWedge X-small, lilac, 36 pcs
602679	FenderWedge & FenderMate Intro-Kit*, 8 & 20 pcs

\*The Intro-Kit contains: FenderWedge x-small 3 pcs, small 3 pcs, medium 1 pc, large 1 pc. FenderMate, left narrow 8 pcs, left regular 2 pcs, right narrow 8 pcs, right regular 2 pcs.

## Confident preparations

High quality preparations made faster without tissue damage.

## Protects the adjacent tooth

The protective shield and wedge helps avoid damage to adjacent tooth.

## Pre-Separation

Facilitates the easy placement of a matrix and the establishment of a good contact.

## Protects gingiva


The protective shield and wedge helps to avoid bleeding and damage to gingival tissue.

# DIRECTA FenderMate®

## FenderMate® Prime

### Fast and safe restorative dentistry with children in mind

FenderMate Prime is a unique combination of a preparation shield and matrix for Class II fillings of primary teeth the "boat-shaped" tip allows the gingiva to be compressed without damage and reduces the risk for bleeding. Neon green FenderMate Prime Long is long & bendable. Neon yellow FenderMate Prime Short is shorter and a more rigid plate for very tight interproximal spaces.

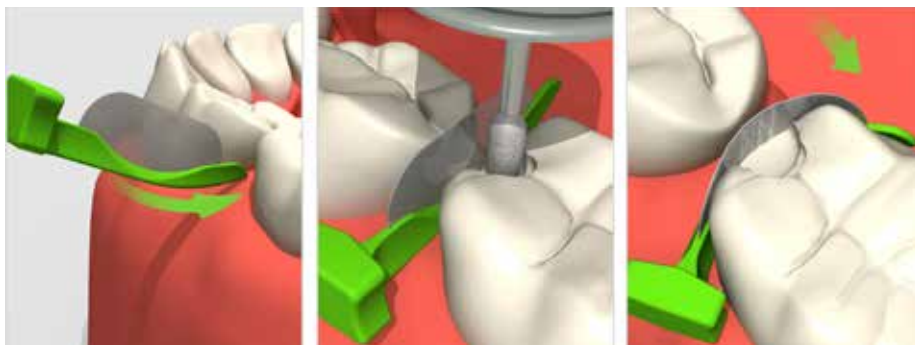

Prod No	Product Description
602815	FenderMate Prime Assortment, 2x18 pcs
602816	FenderMate Prime Short, neon yellow, 18 pcs
602817	FenderMate Prime Long, neon green, 18 pcs

## Designed for primary teeth

The optimized proportions ensure predictable results.

## Easy placement

Facilitates high quality, fast and secure treatments.

## Bendable

Adapts to the curvature of the tooth for a more natural anatomic shape. Minimal finishing required.

## One piece design

Designed to be used for protection during preparation and as a matrix for Class II fillings on primary teeth.



For more information visit [www.directadental.com](http://www.directadental.com) or scan the QR-code with your mobile phone